


PE-Cy5 标记抗人 CD15

目录号	规格
FHC015-02-025	25 次
FHC015-02-100	100 次

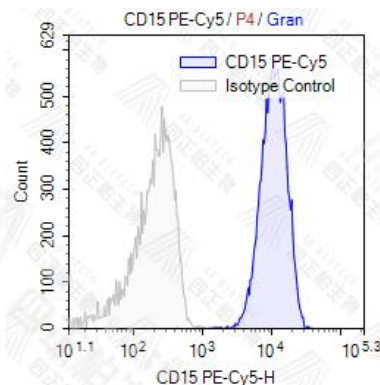

样本：人外周血粒细胞
抗体：PE-Cy5 anti-human CD15(clone 4AW6D3A12) (蓝色)
同型对照：PE-Cy5 mouse IgG1(灰色)

克隆号：	4AW6D3A12
亚型：	mouse IgG1
别名：	Lewis X, 3-FAL, 3-FL, LNFP III, Lex , SSEA-1, X-hapten
缓冲液：	0.2%BSA 和 0.1%Proclin 300 磷酸盐缓冲溶液，PH7.2-7.4
应用：	流式细胞术
推荐用量：	5ul/T
激光器：	蓝色激光 488 nm,绿色激光 532 nm,黄/绿色激光 561 nm
最大吸收/发射波长：	565nm/667nm
反应种属：	人
存储条件：	4 度，避光保存

简介:

CD15 是 3-岩藻糖基-N-乙酰乳糖胺（3-FAL），也被称为 Lewis X、3-FAL、X-hapten 和 SSEA-1。CD15 在粒细胞和单核细胞上表达，也已被证明在朗格汉斯细胞和一些恶性细胞中表达。CD15 已被认为参与粘附、趋化作用、吞噬作用和细菌杀灭活性等过程。

实验方法:

流式细胞术表面染色方案

相关产品:

目录号	名称
FXP055	Sample preservation buffer
FXP001	RBC Lysing Buffer (10x)
FXP005	Cell Staining Buffer
FXP064	Human Fc Receptor Blocking Solution
FMCC001	PE-Cy5 mouse IgG1, Isotype Control

相关文献:

- 1.Wen J, et al. 2020. J Pathol Clin Res. 0.36875.
- 2.Yoshino N, et al. 2000. Exp. Anim. (Tokyo) 49:97. (FC)
- 3.Sander J et al. 2017. Immunity. 47(6):1051-1066 .
- 4.Luoma AM, et al. 2020. Cell. 182(3):655-671.e22.
- 5.Dinh HQ, et al. 2020. Immunity. 53(2):319-334.e6.
- 6.Mathewson ND, et al. 2021. Cell. 184(5):1281-1298.e26.
- 7.Balogh P, et al. 2019. Haematologica. 105:00:00.
- 8.Catapano M, et al. 2020. J Invest Dermatol. 140:816.
- 9.Gong Y, et al. 2020. J Neuroinflammation. 0.845833333.
- 10.Le J, et al. 2020. Immunity. 52(6):1105-1118.e9.

4A Biotech Co., Ltd.

Add: No.Room 201, building 17, No. 8, Jinfeng Road, science and Technology City, Gaoxin District, Suzhou
www.4abio.com , Technical support: tech@4abio.com