

Human GM-CSF ELISA KIT

Catalog Number

CHE0026

Size

96 Tests



Human GM-CSF ELISA KIT

For the quantitative determination of human granulocyte-macrophage colony stimulating factor (GM-CSF) concentrations in cell culture supernates, serum, and plasma. This package insert must be read in its entirety before using this product. If you have questions or experience problems with this product, please contact our Technical Support staff. Our scientists commit themselves to providing rapid and effective help.

**FOR RESEARCH USE ONLY
NOT FOR USE IN DIAGNOSTIC PROCEDURES**

INTRODUCTION

Granulocyte macrophage colony stimulating factor (GM-CSF) is a pleiotropic cytokine with multiple effects on hematopoietic cells. It mobilizes CD34+ progenitor cells into the periphery and stimulates their proliferation, survival and differentiation into neutrophils, monocytes/ macrophages, eosinophils, and myeloid dendritic cells. On these terminally differentiated myeloid cells, GM-CSF is also needed for inducing their effector functions. In addition, GM-CSF has been shown to stimulate the proliferation and differentiation of the erythroid and megakaryocyte progenitor cells .

GM-CSF is produced by a number of different cell types, including keratinocytes, mature and immature NK cells, type II alveolar cells), endothelial cells, monocytes, bone-marrow mesenchymal stem cells, CD4+ and CD8+ T cells, megakaryocytes, B cells, eosinophils, chondrocytes and fibroblasts.

Human GM-CSF cDNA encodes a 144 amino acid (aa) residue precursor protein with a 17 aa putative signal peptide and a 127 aa mature protein. Natural GM-CSF is a monomer that contains both N- and O-linked glycosylation. Mature human GM-CSF shares approximately 55%, 63% and 68% aa sequence homology with mouse, rat and canine. GM-CSF, respectively. Human GM-CSF is not biologically active on mouse cells, but was reported to have some activity on canine cells.

GM-CSF exerts its activity through binding to a high affinity receptor complex consisting of two membrane glycoproteins. The presence of the spliced variants in the heteromeric receptor complex can regulate the functions of the complex.

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for GM-CSF has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any GM-CSF present is bound by the immobilized antibody. Following incubation unbound samples are removed during a wash step, and then a detection antibody specific for GM-CSF is added to the wells and binds to the combination of capture antibody- GM-CSF in sample. Following a wash to remove any unbound combination, and enzyme conjugate is added to the wells. Following incubation and wash steps a substrate is added. A coloured product is formed in proportion to the amount of GM-CSF present in the sample. The reaction is terminated by

addition of acid and absorbance is measured at 450nm. A standard curve is prepared from seven GM-CSF standard dilutions and GM-CSF sample concentration determined.

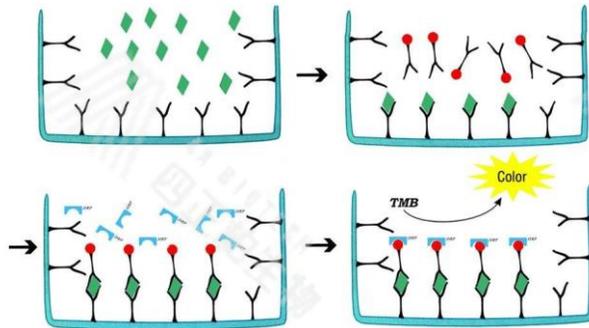


Figure 1: Schematic diagram of the assay

REAGENTS

1. Aluminium pouches with a Microwell Plate coated with monoclonal antibody to human GM-CSF (8×12)
2. 2 vials human GM-CSF Standard lyophilized, 2000 pg/ml upon reconstitution
3. 2 vials concentrated Biotin-Conjugate anti-human GM-CSF monoclonal antibody
4. 2 vials Streptavidin-HRP solution,
5. 1 bottle Standard /sample Diluent
6. 1 bottle Biotin-Conjugate antibody Diluent
7. 1 bottle Streptavidin-HRP Diluent
8. 1 bottle Wash Buffer Concentrate 20x (PBS with 1% Tween-20)
9. 1 vial Substrate Solution
10. 1 vial Stop Solution
11. 4 pieces Adhesive Films
12. package insert

NOTE: [96 Tests]

STORAGE

Table 1: Storage of the kit

| | | |
|---|---|---|
| Unopened Kit | Store at 2 - 8° C. Do not use past kit expiration date. | |
| | Standard /sample Diluent | May be stored for up to 1 month at 2 - 8° C.** |
| | Concentrated Biotin-Conjugate | |
| | Streptavidin-HRP solution | |
| | Biotin-Conjugate antibody Diluent | |
| | Streptavidin-HRP Diluent | |
| | Wash Buffer Concentrate 20x | |
| | Substrate Solution | |
| | Stop Solution | |
| Opened/ Reconstituted Reagents | Standard | Aliquot and store for up to 1 month at ≤-20°C. Avoid repeated freeze-thaw cycles. Diluted standard shall not be reused. |
| | Microplate Wells | Return unused wells to the foil pouch containing the desiccant pack, reseal along entire edge of zip-seal. May be stored for up to 1 month at 2 - 8° C.** |

**Provided this is within the expiration date of the kit.

THE REQUIRED ITEMS (not provided, but can help to buy):

1. Microplate reader (450nm).
2. Micro-pipette and tips: 0.5-10, 2-20, 20-200, 200-1000ul.
3. 37 °C incubator, double-distilled water or deionized water, coordinate paper, graduated cylinder.

PRECAUTIONS FOR USE

1. Store kit reagents between 2°C and 8°C. After use all reagents should be immediately returned to cold storage (2°C to 8°C).
2. Please perform simple centrifugation to collect the liquid before use.

3. To avoid cross contamination, please use disposable pipette tips.
4. The Stop Solution suggested for use with this kit is an acid solution. Wear eye, hand, face, and clothing protection when using this material. Avoid contact of skin or mucous membranes with kit reagents or specimens. In the case of contact with skin or eyes wash immediately with water.
5. Use clean, dedicated reagent trays for dispensing the washing liquid, conjugate and substrate reagent. Mix all reagents and samples well before use.
6. After washing microtiter plate should be fully pat dried. Do not use absorbent paper directly into the enzyme reaction wells.
7. Do not mix or substitute reagents with those from other lots or other sources. Do not use kit reagents beyond expiration date on label.
8. Each sample, standard, blank and optional control samples should be assayed in duplicate or triplicate.
9. Adequate mixing is very important for good result. Use a mini-vortexer at the lowest frequency or Shake by hand at 10min interval when there is no vortexer.
10. Avoid microtiter plates drying during the operation.
11. Dilute samples at the appropriate multiple, and make the sample values fall within the standard curve. If samples generate values higher than the highest standard, dilute the samples and repeat the assay.
12. Any variation in standard diluent, operator, pipetting technique, washing technique, incubation time and temperature, and kit age can cause variation in binding.
13. This method can effectively eliminate the interference of the soluble receptors, binding proteins and other factors in biological samples.

SAMPLE COLLECTION AND STORAGE

1. **Cell Culture Supernates** - Remove particulates by centrifugation.
2. **Serum** - Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 15 minutes at approximately 1000 x g. Remove serum, avoid hemolysis and high blood lipid samples.
3. **Plasma** - Recommended EDTA as an anticoagulant in plasma. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection.

4. Assay immediately or aliquot and store samples at -20°C. Avoid repeated freeze-thaw cycles.
5. Dilute samples at the appropriate multiple (recommended to do pre-test to determine the dilution factor).

Note: The normal human serum or plasma samples are suggested to make a 1:2 dilution.

REAGENT PREPARATION

1. Bring all reagents to room temperature before use.
2. **Wash Buffer** - Dilute 10mL of Wash Buffer Concentrate into deionized or distilled water to prepare 200mL of Wash Buffer. If crystals have formed in the concentrate Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.
3. **Standard** - Reconstitute the Standard with 1.0mL of Standard /sample Diluent. This reconstitution produces a stock solution of 2000 pg /mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions.

Pipette 500 μ L of Standard/sample Diluent into the 1000 pg/mL tube and the remaining tubes. Use the stock solution to produce a 2-fold dilution series (below). Mix each tube thoroughly and change pipette tips between each transfer. The 1000 pg/mL standard serves as the high standard. The Standard/ sample Diluent serves as the zero standard (0 pg/mL).

If you do not run out of re-melting standard, store it at -20°C. Diluted standard shall not be reused.

4. Working solution of Biotin-Conjugate anti-human GM-CSF monoclonal antibody:
Make a 1:100 dilution of the concentrated Biotin-Conjugate solution with the Biotin-Conjugate antibody Diluent in a clean plastic tube.

The working solution should be used within one day after dilution.

5. Working solution of Streptavidin-HRP: Make a 1:100 dilution of the concentrated Streptavidin-HRP solution with the Streptavidin-HRP Diluent in a clean plastic tube.

The working solution should be used within one day after dilution.

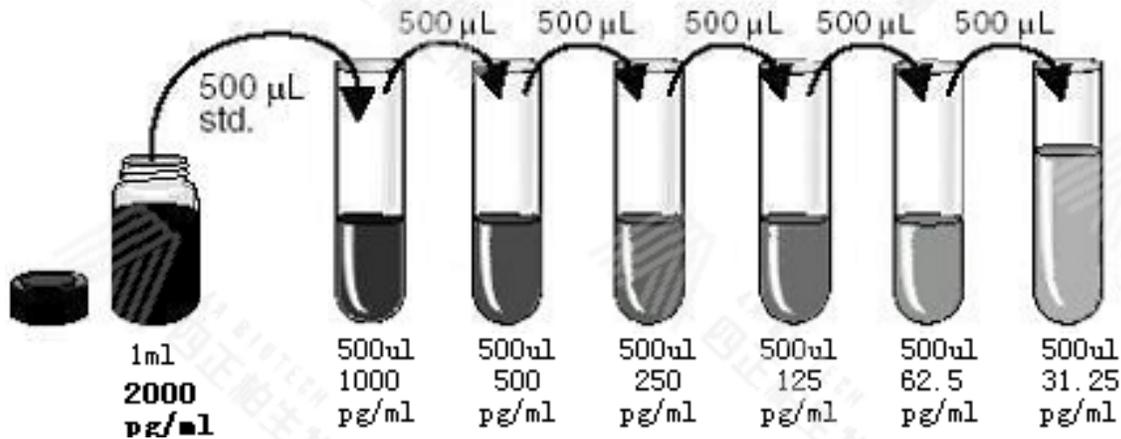


Figure 2: Preparation of GM-CSF standard dilutions

GENERAL ELISA PROTOCOL

1. Prepare all reagents and working standards as directed in the previous sections.
2. Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running blanks and standards. Remove extra microwell strips from holder and store in foil bag with the desiccant provided at 2-8°C sealed tightly.
3. Add 100µL of Standard, control, or sample, per well. Cover with the adhesive strip provided. Incubate for 1.5 hours at 37°C.
4. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with Wash Buffer (350µL) using a squirt bottle, manifold dispenser or auto-washer. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
5. Add 100 µL of the working solution of Biotin-Conjugate to each well. Cover with a new adhesive strip and incubate 1 hours at 37°C.
6. Repeat the aspiration/wash as in step 3.
7. Add 100 µL of the working solution of Streptavidin-HRP to each well. Cover with a new adhesive strip and incubate for 30 minutes at 37°C. Avoid placing the plate in direct light.
8. Repeat the aspiration/wash as in step 3.

9. Add 100 μ L of Substrate Solution to each well. Incubate for 10-20 minutes at 37°C. Avoid placing the plate in direct light.
10. Add 100 μ L of Stop Solution to each well. Gently tap the plate to ensure thorough mixing.
11. Determine the optical density of each well immediately, using a microplate reader set to 450 nm.(optionally 630nm as the reference wave length;610-650nm is acceptable)

ASSAY PROCEDURE SUMMARY

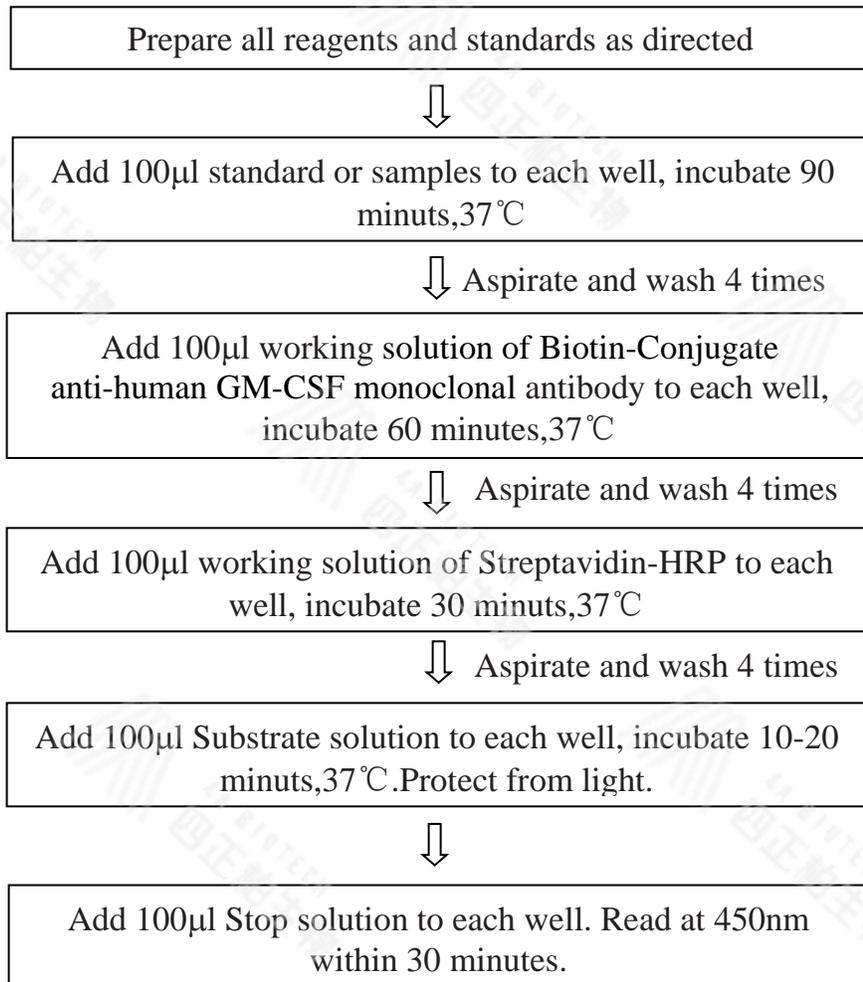


Figure 3: Assay procedure summary

TECHNICAL HINTS

1. When mixing or reconstituting protein solutions, always avoid foaming.
2. To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
3. To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary.
4. Substrate Solution should remain colorless until added to the plate. Stop Solution should be added to the plate in the same order as the Substrate Solution. Keep Substrate Solution protected from light. Substrate Solution should change from colorless to gradations of blue.
5. A standard curve should be generated for each set of samples assayed. According to the content of tested factors in the sample, appropriate diluted or concentrated samples, it is best to do pre-experiment.

CALCULATION OF RESULTS

1. Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density.
2. Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph.
3. The data may be linearized by plotting the log of the GM-CSF concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.
4. This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Table 2: Typical data using the GM-CSF ELISA (Measuring wavelength: 450nm, Reference wavelength: 630nm)

| Standardized (pg/ml) | OD. | OD. | Average | Corrected |
|----------------------|-------|-------|---------|-----------|
| 0 | 0.089 | 0.092 | 0.091 | — |
| 15.625 | 0.143 | 0.146 | 0.145 | 0.054 |
| 31.25 | 0.238 | 0.224 | 0.231 | 0.140 |
| 62.5 | 0.343 | 0.346 | 0.345 | 0.254 |
| 125 | 0.582 | 0.591 | 0.586 | 0.495 |
| 250 | 0.927 | 0.914 | 0.921 | 0.830 |
| 500 | 1.542 | 1.553 | 1.548 | 1.457 |
| 1000 | 2.318 | 2.310 | 2.314 | 2.223 |

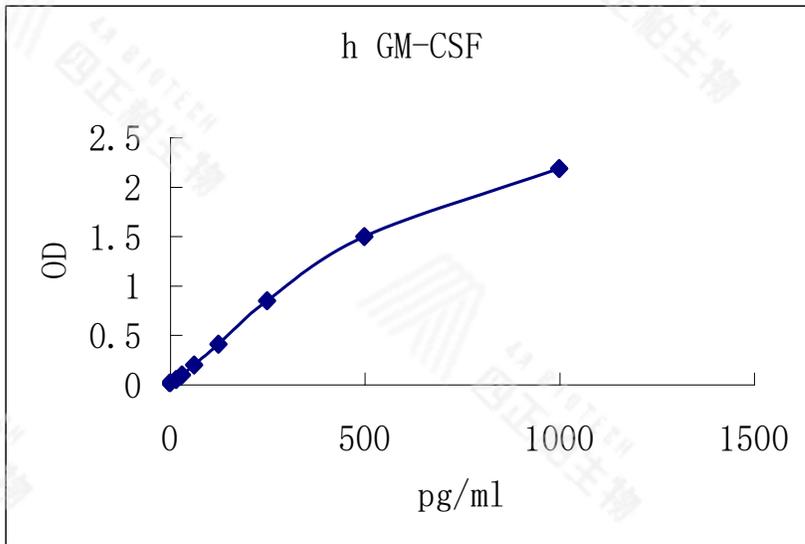


Figure 4: Representative standard curve for GM-CSF ELISA. GM-CSF was diluted in serial two-fold steps in Sample Diluent.

Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

SENSITIVITY, SPECIFICITY AND REPEATABILITY

1. **REPEATABILITY:** The coefficient of variation of both intra-assay and inter-assay were less than 10%.
2. **SENSITIVITY:** The minimum detectable dose was 4pg/mL.
3. **SPECIFICITY:** This assay recognizes both natural and recombinant human GM-CSF. The factors listed below were prepared at 10ng/ml in Standard /sample Diluent and assayed for cross-reactivity and no significant cross-reactivity or interference was observed.

Table 3: Factors assayed for cross-reactivity

| Recombinant human | Recombinant mouse | Other proteins |
|-------------------|-------------------|------------------------|
| G-CSF | GM-CSF | bovine FGF acidic |
| M-CSF | IL-1 α | bovine FGF basic |
| FGF-4 | IL-1 β | human PDGF |
| FGF-5 | IL-3 | porcine PDGF |
| FGF-6 | IL-4 | Human TGF- β 1 |
| IL-1 α | IL-5 | porcine TGF- β 1 |
| IL-1 β | IL-6 | |
| IL-3 | IL-7 | |
| IL-5 | IL-9 | |
| IL-6 | IL-10 | |
| IL-7 | TNF- α | |
| IL-8 | | |
| IL-10 | | |
| IL-12 | | |
| VEGF | | |

REFERENCES

1. Cantrell, M.A. et al. (1985) Proc. Natl. Acad. Sci. USA 82:6250.
2. Wong, G.C. et al. (1985) Science 228:810.
3. Lee, F. et al. (1985) Proc. Natl. Acad. Sci. USA 82:4360.
4. Smith, L.R. et al. (1994) Immunogenetics 39:80.
5. Gough, N.M. et al. (1985) EMBO J. 4:645.
6. Nash, R.A. et al. (1991) Blood 78:930.

7. Guthridge, M.A. et al. (1998) Stem Cells 16:301.
8. Gearing, D.P. et al. (1989) EMBO J. 8:3667.
9. Hayashida, K. et al. (1990) Proc. Natl. Acad. Sci. USA 87:9655.
10. Till, S. et al. (1995) Eur. J. Immunol. 25:2727.
11. Harris, R.J. et al. (2000) J. Immunol. 164:3887.
12. Loza, M.J. et al. (2002) Blood 99:1273.

RELATED PRODUCTS

Table 4: Related products

| Products name | Catalog number | size |
|--------------------------------|----------------|---------|
| Human IL-1 α ELISA Kit | CHE0101 | 48T/96T |
| Human IL-1 β ELISA Kit | CHE0001 | 48T/96T |
| Human IL-2 ELISA Kit | CHE0003 | 48T/96T |
| Human IL-3 ELISA Kit | CHE0028 | 48T/96T |
| Human IL-4 ELISA Kit | CHE0005 | 48T/96T |
| Human IL-6 ELISA Kit | CHE0009 | 48T/96T |
| Human IL-8 ELISA Kit | CHE0011 | 48T/96T |
| Human IL-10 ELISA Kit | CHE0013 | 48T/96T |
| Human IL-12p40 ELISA Kit | CHE0016 | 48T/96T |
| Human IL-12p70 ELISA Kit | CHE0015 | 48T/96T |
| Human IFN- γ ELISA Kit | CHE0017 | 48T/96T |
| Human TNF- α ELISA Kit | CHE0019 | 48T/96T |
| Human APO-1/FAS ELISA Kit | CHE0024 | 48T/96T |
| Human TGF- β 1 ELISA Kit | CHE0029 | 48T/96T |
| Human MCP-1 ELISA Kit | CHE0103 | 48T/96T |
| Human EGF ELISA Kit | CHE0044 | 48T/96T |
| Human G-CSF ELISA Kit | CHE0025 | 48T/96T |
| Human GM-CSF ELISA Kit | CHE0026 | 48T/96T |
| Human sICAM-1 ELISA Kit | CHE0052 | 48T/96T |
| Human Leptin ELISA Kit | CHE0053 | 48T/96T |
| Human VEGF ELISA Kit | CHE0043 | 48T/96T |

If you have any questions, please tell us!