

Human sICAM-1 ELISA KIT

Catalog Number

CHE0052

Size

96 Tests



Human sICAM-1 ELISA KIT

For the quantitative determination of human soluble Intercellular Adhesion Molecule 1 (sICAM-1) concentrations in cell culture supernates, serum, and plasma. This package insert must be read in its entirety before using this product. If you have questions or experience problems with this product, please contact our Technical Support staff. Our scientists commit themselves to providing rapid and effective help.

**FOR RESEARCH USE ONLY
NOT FOR USE IN DIAGNOSTIC PROCEDURES**

INTRODUCTION

Intercellular Adhesion Molecule 1 (ICAM-1), also known as CD54, is a nearly ubiquitous transmembrane glycoprotein that plays a key role in leukocyte migration and activation. Human ICAM-1 contains five Ig-like domains in its extracellular domain (ECD) and associates into non-covalently linked dimers. Soluble forms of monomeric and dimeric ICAM-1 (sICAM-1) can be generated via proteolytic cleavage by cathepsin G, elastase, MMP-9, MMP-14/MT1-MMP, and TACE/ADAM17. In the mouse, alternate splicing generates isoforms that lack particular Ig-like domains and are differentially sensitive to proteolysis. Within the ECD, human ICAM-1 shares 53% amino acid sequence identity with mouse and rat ICAM-1.

The principal binding partners of ICAM-1 are the leukocyte integrins LFA-1 (CD11a/CD18) and Mac-1 (CD11b/CD18). The multivalency of dimeric ICAM-1 increases its strength of interaction with LFA-1. ICAM-1 also binds several non-integrin ligands including CD43/sialoporphin, fibrinogen, hyaluronan, rhinoviruses, and Plasmodium falciparum-infected erythrocytes. At sites of inflammation, ICAM-1 is upregulated on endothelial and epithelial cells where it mediates the adhesion and paracellular migration of leukocytes expressing activated LFA-1 and Mac-1. ICAM-1 ligation prolongs antigen presentation by dendritic cells and promotes T cell proliferation and cytokine release. ICAM-1 activation also participates in angiogenesis, wound healing, and bone metabolism.

Soluble ICAM-1 has been reported in serum, cerebrospinal fluid, urine, and bronchoalveolar lavage fluid. Elevated levels of sICAM-1 in these fluids are associated with cardiovascular disease, type 2 diabetes, organ transplant dysfunction, oxidant stress, abdominal fat mass, hypertension, liver disease, and certain malignancies. sICAM-1 promotes angiogenesis and serves as an indicator of vascular endothelial cell activation or damage. It also functions as an inhibitor of transmembrane ICAM-1 mediated activities such as monocyte adhesion to activated endothelial cells and sensitivity of tumor cells to NK cell-mediated lysis.

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for sICAM-1 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any sICAM-1 present is bound by the immobilized antibody. Following incubation unbound samples are removed during a wash step, and then a detection antibody specific for sICAM-1 is added to the wells and binds to the combination of capture antibody- sICAM-1 in sample.

Following a wash to remove any unbound combination, and enzyme conjugate is added to the wells. Following incubation and wash steps a substrate is added. A coloured product is formed in proportion to the amount of sICAM-1 present in the sample. The reaction is terminated by addition of acid and absorbance is measured at 450nm. A standard curve is prepared from seven sICAM-1 standard dilutions and sICAM-1 sample concentration determined.

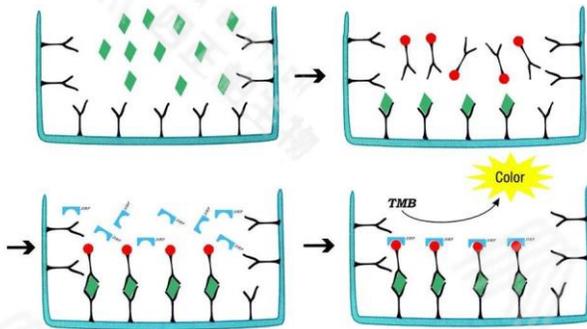


Figure 1: Schematic diagram of the assay

REAGENTS

1. Aluminium pouches with a Microwell Plate coated with antibody to human sICAM-1 (8×12)
2. 2 vials human sICAM-1 Standard lyophilized, 4000 pg/ml upon reconstitution
3. 2 vials concentrated Biotin-Conjugate anti-human sICAM-1 antibody
4. 2 vials Streptavidin-HRP solution,
5. 4 bottle Standard /sample Diluent
6. 1 bottle Biotin-Conjugate antibody Diluent
7. 1 bottle Streptavidin-HRP Diluent
8. 1 bottle Wash Buffer Concentrate 20x (PBS with 1% Tween-20)
9. 1 vial Substrate Solution
10. 1 vial Stop Solution
11. 4 pieces Adhesive Films
12. package insert

NOTE: [96 Tests]

STORAGE

Table 1: Storage of the kit

Unopened Kit	Store at 2 - 8° C. Do not use past kit expiration date.	
Opened/ Reconstituted Reagents	Standard /sample Diluent	May be stored for up to 1 month at 2 - 8° C.**
	Concentrated Biotin-Conjugate	
	Streptavidin-HRP solution	
	Biotin-Conjugate antibody Diluent	
	Streptavidin-HRP Diluent	
	Wash Buffer Concentrate 20x	
	Substrate Solution	
	Stop Solution	
	Standard	Aliquot and store for up to 1 month at ≤-20°C. Avoid repeated freeze-thaw cycles. Diluted standard shall not be reused.
	Microplate Wells	Return unused wells to the foil pouch containing the desiccant pack, reseal along entire edge of zip-seal. May be stored for up to 1 month at 2 - 8° C.**

**Provided this is within the expiration date of the kit.

THE REQUIRED ITEMS (not provided, but can help to buy):

1. Microplate reader (450nm).
2. Micro-pipette and tips: 0.5-10, 2-20, 20-200, 200-1000ul.
3. 37 °C incubator, double-distilled water or deionized water, coordinate paper, graduated cylinder.

PRECAUTIONS FOR USE

1. Store kit reagents between 2°C and 8°C. After use all reagents should be immediately returned to cold storage (2°C to 8°C).
2. Please perform simple centrifugation to collect the liquid before use.

3. To avoid cross contamination, please use disposable pipette tips.
4. The Stop Solution suggested for use with this kit is an acid solution. Wear eye, hand, face, and clothing protection when using this material. Avoid contact of skin or mucous membranes with kit reagents or specimens. In the case of contact with skin or eyes wash immediately with water.
5. Use clean, dedicated reagent trays for dispensing the washing liquid, conjugate and substrate reagent. Mix all reagents and samples well before use.
6. After washing microtiter plate should be fully pat dried. Do not use absorbent paper directly into the enzyme reaction wells.
7. Do not mix or substitute reagents with those from other lots or other sources. Do not use kit reagents beyond expiration date on label.
8. Each sample, standard, blank and optional control samples should be assayed in duplicate or triplicate.
9. Adequate mixing is very important for good result. Use a mini-vortexer at the lowest frequency or Shake by hand at 10min interval when there is no vortexer.
10. Avoid microtiter plates drying during the operation.
11. Dilute samples at the appropriate multiple, and make the sample values fall within the standard curve. If samples generate values higher than the highest standard, dilute the samples and repeat the assay.
12. Any variation in standard diluent, operator, pipetting technique, washing technique, incubation time and temperature, and kit age can cause variation in binding.
13. This method can effectively eliminate the interference of the soluble receptors, binding proteins and other factors in biological samples.

SAMPLE COLLECTION AND STORAGE

1. **Cell Culture Supernates** - Remove particulates by centrifugation.
2. **Serum** - Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 15 minutes at approximately 1000 x g. Remove serum, avoid hemolysis and high blood lipid samples.
3. **Plasma** - Recommended EDTA as an anticoagulant in plasma. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection.

4. Assay immediately or aliquot and store samples at -20°C. Avoid repeated freeze-thaw cycles.
5. Dilute samples at the appropriate multiple (recommended to do pre-test to determine the dilution factor).

Note: The normal human serum or plasma samples are suggested to make a 1:150-1:200 dilution.

REAGENT PREPARATION

1. Bring all reagents to room temperature before use.
2. **Wash Buffer** - Dilute 10mL of Wash Buffer Concentrate into deionized or distilled water to prepare 200mL of Wash Buffer. If crystals have formed in the concentrate Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.
3. **Standard** - Reconstitute the Standard with 1.0mL of Standard /sample Diluent. This reconstitution produces a stock solution of 4000 pg /mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions.

Pipette 500 μ L of Standard/sample Diluent into the 2000 pg/mL tube and the remaining tubes. Use the stock solution to produce a 2-fold dilution series (below). Mix each tube thoroughly and change pipette tips between each transfer. The 4000 pg/mL standard serves as the high standard. The Standard/ sample Diluent serves as the zero standard (0 pg/mL).

If you do not run out of re-melting standard, store it at -20°C. Diluted standard shall not be reused.

4. Working solution of Biotin-Conjugate anti-human sICAM-1 antibody: Make a 1:100 dilution of the concentrated Biotin-Conjugate solution with the Biotin-Conjugate antibody Diluent in a clean plastic tube.

The working solution should be used within one day after dilution.

5. Working solution of Streptavidin-HRP: Make a 1:100 dilution of the concentrated Streptavidin-HRP solution with the Streptavidin-HRP Diluent in a clean plastic tube.

The working solution should be used within one day after dilution.

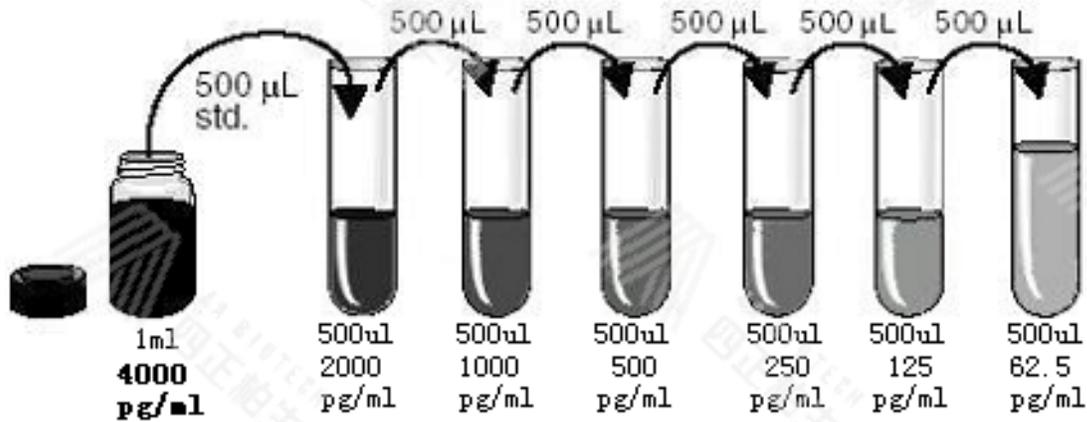


Figure 2: Preparation of sICAM-1 standard dilutions

GENERAL ELISA PROTOCOL

1. Prepare all reagents and working standards as directed in the previous sections.
2. Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running blanks and standards. Remove extra microwell strips from holder and store in foil bag with the desiccant provided at 2-8°C sealed tightly.
3. Add 100 µL of Standard, control, or sample, per well. Cover with the adhesive strip provided. Incubate for 1.5 hours at 37°C.
4. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with Wash Buffer (350 µL) using a squirt bottle, manifold dispenser or auto-washer. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
5. Add 100 µL of the working solution of Biotin-Conjugate to each well. Cover with a new adhesive strip and incubate 1 hours at 37°C.
6. Repeat the aspiration/wash as in step 4.
7. Add 100 µL of the working solution of Streptavidin-HRP to each well. Cover with a new adhesive strip and incubate for 30 minutes at 37°C. Avoid placing the plate in direct light.
8. Repeat the aspiration/wash as in step 4.
9. Add 100 µL of Substrate Solution to each well. Incubate for 10-20

minutes at 37°C. Avoid placing the plate in direct light.

10. Add 100 μ L of Stop Solution to each well. Gently tap the plate to ensure thorough mixing.
11. Determine the optical density of each well immediately, using a microplate reader set to 450 nm.(optionally 630nm as the reference wave length;610-650nm is acceptable)

ASSAY PROCEDURE SUMMARY

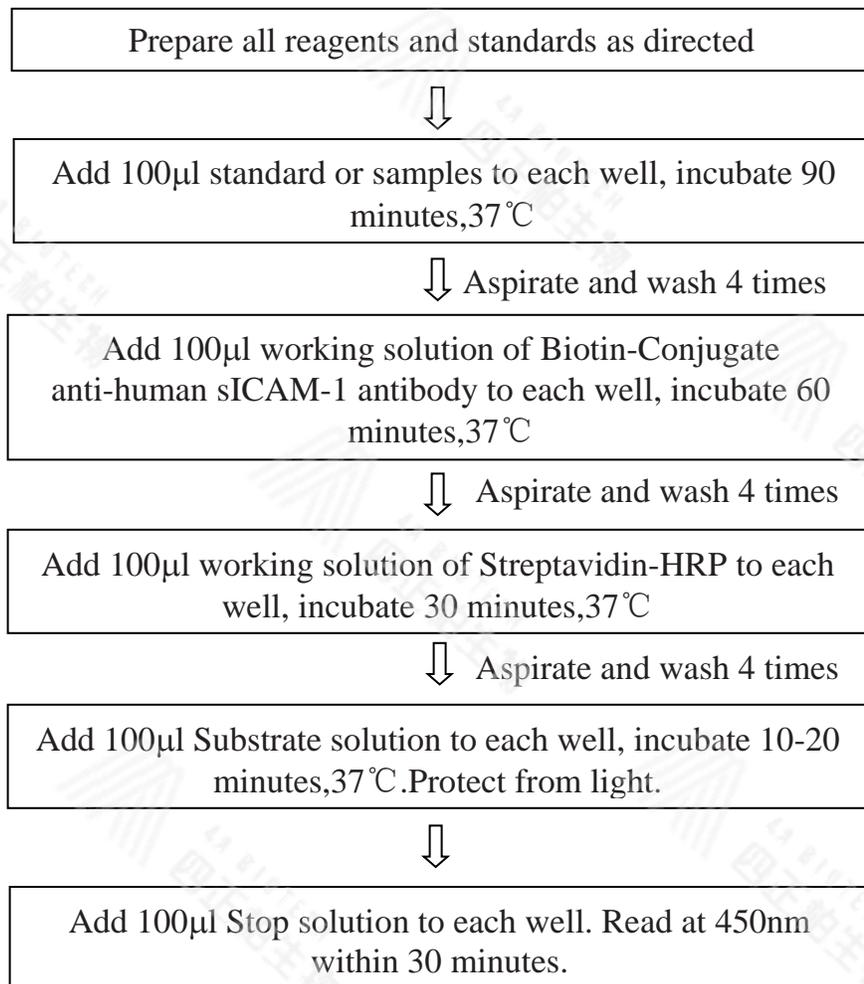


Figure 3: Assay procedure summary

TECHNICAL HINTS

1. When mixing or reconstituting protein solutions, always avoid foaming.
2. To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
3. To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary.
4. Substrate Solution should remain colorless until added to the plate. Stop Solution should be added to the plate in the same order as the Substrate Solution. Keep Substrate Solution protected from light. Substrate Solution should change from colorless to gradations of blue.
5. A standard curve should be generated for each set of samples assayed. According to the content of tested factors in the sample, appropriate diluted or concentrated samples, it is best to do pre-experiment.

CALCULATION OF RESULTS

1. Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density.
2. Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph.
3. The data may be linearized by plotting the log of the sICAM-1 concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.
4. This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Table 2: Typical data using the sICAM-1 ELISA (Measuring wavelength:450nm,Reference wavelength:630nm)

Standard (pg/ml)	OD.	OD.	Average	Corrected
0	0.043	0.042	0.042	—
62.5	0.146	0.142	0.144	0.102
125	0.213	0.210	0.211	0.169
250	0.349	0.345	0.347	0.305
500	0.599	0.597	0.598	0.556
1000	0.949	0.946	0.947	0.905
2000	1.557	1.555	1.556	1.514
4000	1.947	1.946	1.946	1.904

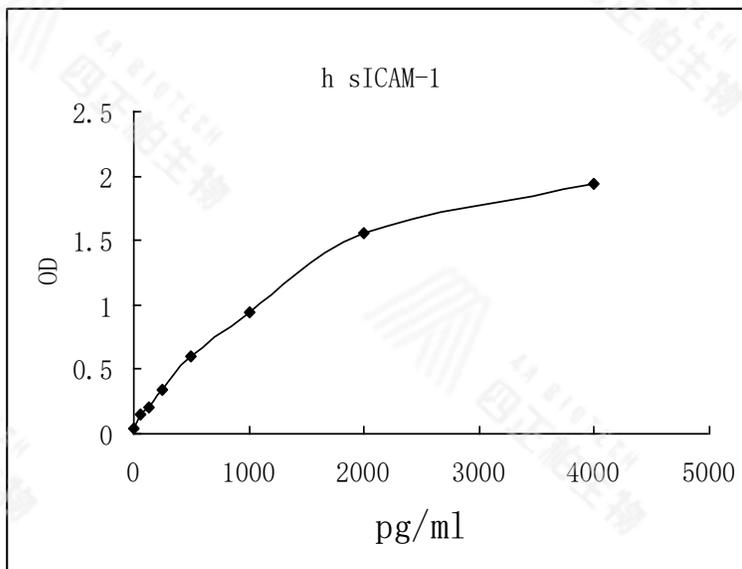


Figure 4: Representative standard curve for sICAM-1 ELISA. sICAM-1 was diluted in serial two-fold steps in Sample Diluent.

Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

SENSITIVITY, SPECIFICITY AND REPEATABILITY

1. **REPEATABILITY:** The coefficient of variation of both intra-assay and inter-assay were less than 10%.
2. **SENSITIVITY:** The minimum detectable dose was 31pg/mL.
3. **SPECIFICITY:** This assay recognizes both natural and recombinant human sICAM-1. The factors listed below were prepared at 50ng/ml in Standard /sample Diluent and assayed for cross-reactivity and no significant cross-reactivity or interference was observed.

Table 3: Factors assayed for cross-reactivity

Recombinant human	Recombinant mouse
VCAM-1	JAM-A
L-Selectin	JAM-B
Siglec-2	JAM-C
Siglec-3	ICAM-1
Siglec-5	ICAM-2
ICAM-2	ICAM-5
ICAM-3	CHL-1
ICAM-5	VCAM-1
JAM-A	P-Selectin
JAM-B	L-Selectin
JAM-C	
CHL-1	
BCAM	
CD31	

REFERENCES

1. Hogg, N. et al. (1991) Chem. Immunol. 50:98.
2. Witkowska, A.M. and M.H. Borawska (2004) Eur. Cytokine Netw. 15:91.
3. Staunton, D.E. et al. (1988) Cell 52:925.
4. Reilly, P.L. et al. (1995) J. Immunol. 155:529.
5. Robledo, O. et al. (2003) Eur. J. Immunol. 33:1351.
6. Sithu, S.D. et al. (2007) J. Biol. Chem. 282:25010.
7. Fiore, E. et al. (2002) Oncogene 21:5213.
8. Tsakadze, N.L. et al. (2006) J. Biol. Chem. 281:3157.
9. Jun, C.D. et al. (2001) Proc. Natl. Acad. Sci. 98:6830.

RELATED PRODUCTS

Table 4: Related products

Products name	Catalog number	size
Human IL-1 α ELISA Kit	CHE0101	48T/96T
Human IL-1 β ELISA Kit	CHE0001	48T/96T
Human IL-2 ELISA Kit	CHE0003	48T/96T
Human IL-3 ELISA Kit	CHE0028	48T/96T
Human IL-4 ELISA Kit	CHE0005	48T/96T
Human IL-6 ELISA Kit	CHE0009	48T/96T
Human IL-8 ELISA Kit	CHE0011	48T/96T
Human IL-10 ELISA Kit	CHE0013	48T/96T
Human IL-12p40 ELISA Kit	CHE0016	48T/96T
Human IL-12p70 ELISA Kit	CHE0015	48T/96T
Human IFN- γ ELISA Kit	CHE0017	48T/96T
Human TNF- α ELISA Kit	CHE0019	48T/96T
Human APO-1/FAS ELISA Kit	CHE0024	48T/96T
Human TGF- β 1 ELISA Kit	CHE0029	48T/96T
Human MCP-1 ELISA Kit	CHE0103	48T/96T
Human EGF ELISA Kit	CHE0044	48T/96T
Human G-CSF ELISA Kit	CHE0025	48T/96T
Human GM-CSF ELISA Kit	CHE0026	48T/96T
Human sICAM-1 ELISA Kit	CHE0052	48T/96T
Human Leptin ELISA Kit	CHE0053	48T/96T
Human VEGF ELISA Kit	CHE0043	48T/96T

If you have any questions, please tell us!